

# Product Catalogue

## NAÍDA CI M & SKY CI M

Sound Processors and Accessories



NAÍDA CI M | SKY CI M







# CONNECT WITH THE MOMENTS YOU LOVE

<b>MARVEL CI SOUND PROCESSORS</b> .....	4
Colour Options.....	5
<b>SLIM HP</b> .....	6
Headpieces .....	6
Colour Caps.....	6
Slim HP Magnets .....	7
<b>WATERPROOF ACCESSORIES</b> .....	8
Waterproof Battery & Slim HP AquaMic .....	8
Colour Options.....	9
<b>POWER OPTIONS</b> .....	10
Colour Options.....	11
Charger .....	11
<b>EARHOOK</b> .....	12
<b>ACOUSTIC EARHOOK</b> .....	12
<b>T-MIC MICROPHONE</b> .....	12
<b>WEARING AND RETENTION OPTIONS</b> .....	13
<b>ADDITIONAL ACCESSORIES</b> .....	14
AB Remote app .....	15
<b>CONNECTIVITY</b> .....	16
Roger Accessories .....	16
Phonak Accessories .....	17
<b>BIMODAL SOLUTION</b> .....	18
Phonak Naída Link M .....	20
Phonak Sky Link M .....	20
<b>TECHNICAL OVERVIEW</b> .....	22

# MARVEL CI SOUND PROCESSORS DESIGNED WITH YOU IN MIND

Naída™ CI M and Sky CI™ M sound processors are part of an innovative, advanced cochlear implant system that delivers a brilliant hearing experience. This means powerful sound and excellent hearing in any environment, along with flexibility and durability that allows you or your children to meet all of life's adventures head-on.

Naída CI M and Sky CI M offer a wide **variety of features** designed to improve hearing across all listening situations using **AutoSense OS™ 3.0 operating system and direct streaming** from a multitude of devices. Rest assured—every aspect of the system is designed to be seamlessly operated, providing the ultimate ease of use. Please refer to the Technical Overview (page 22) for additional information on product and feature differentiation.

## COMPONENTS

1. M T-Mic™ Microphone with Cover
2. Sound Processor
3. Front Microphone
4. Multi-function Button with Status LED
5. Back Microphone
6. M Battery
7. Slim HP with Integrated Cable




### NAÍDA CI M90

### SKY CI M90

		NAÍDA CI M90	SKY CI M90
ALPINE WHITE	●	N/A	CI-5295-110
SAND BEIGE	●	CI-5293-120	CI-5295-120
CHESTNUT	●	CI-5293-130	CI-5295-130
SILVER GREY	●	CI-5293-140	CI-5295-140
VELVET BLACK	●	CI-5293-150	CI-5295-150
CARIBBEAN PIRATE	●	N/A	CI-5295-190
PRECIOUS PINK	●	N/A	CI-5295-250
BLUE OCEAN	●	N/A	CI-5295-260
MAJESTY PURPLE	●	N/A	CI-5295-270
LAVA RED	●	N/A	CI-5295-280
RESIN BEIGE*	●	CI-5293-240	CI-5295-240

\* These unpainted sound processors are intended for hypo-allergenic patients who could have allergic reactions to paint. Made to order, please contact your local AB representative prior to ordering.



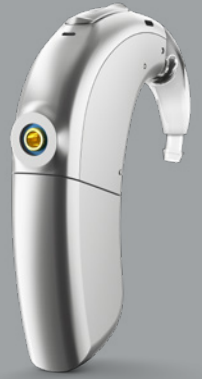














## SLIM HEADPIECE

Naída CI M and Sky CI M feature an all-new, **slim headpiece** which offers comfort and convenience. With a low profile, it stays out of the way so you can enjoy an active lifestyle.



















SLIM HP 9 CM (3.5") BEIGE <sup>1</sup>		CI-5320-001
SLIM HP 9 CM (3.5") BLACK <sup>1</sup>		CI-5320-002
SLIM HP 9 CM (3.5") BROWN <sup>1</sup>		CI-5320-003
SLIM HP 9 CM (3.5") GREY <sup>1</sup>		CI-5320-004
SLIM HP 9 CM (3.5") WHITE <sup>1</sup>		CI-5320-005
SLIM HP 11 CM (4.25") BEIGE <sup>1</sup>		CI-5320-006
SLIM HP 11 CM (4.25") BLACK <sup>1</sup>		CI-5320-007
SLIM HP 11 CM (4.25") BROWN <sup>1</sup>		CI-5320-008
SLIM HP 11 CM (4.25") GREY <sup>1</sup>		CI-5320-009
SLIM HP 11 CM (4.25") WHITE <sup>1</sup>		CI-5320-010
SLIM HP 14 CM (5.5") BEIGE <sup>1</sup>		CI-5320-011
SLIM HP 14 CM (5.5") BLACK <sup>1</sup>		CI-5320-012
SLIM HP M 30.5 CM (12") BLACK <sup>2</sup>		CI-5322-001

1. For on-ear-use (no integrated microphone)

2. For off-ear-use (with integrated microphone) information

## SLIM HP COLOUR

COLOUR CAP LAVA RED		CI-7131-001
COLOUR CAP BLUE OCEAN		CI-7131-002
COLOUR CAP MAJESTY PURPLE		CI-7131-003
COLOUR CAP SAND BEIGE		CI-7131-004
COLOUR CAP CHESTNUT		CI-7131-005
COLOUR CAP SILVER GREY		CI-7131-006
COLOUR CAP VELVET BLACK		CI-7131-007
COLOUR CAP CARIBBEAN PIRATE		CI-7131-008
COLOUR CAP PRECIOUS PINK		CI-7131-009
COLOUR CAP ALPINE WHITE		CI-7131-010
COLOUR CAP RESIN BLACK*		CI-7131-011
COLOUR CAP SAND BEIGE LARGE		CI-7132-001
COLOUR CAP CHESTNUT LARGE		CI-7132-002
COLOUR CAP SILVER GREY LARGE		CI-7132-003
COLOUR CAP VELVET BLACK LARGE		CI-7132-004
COLOUR CAP RESIN BLACK LARGE*		CI-7132-005
SLIM HP COLOUR CAP REMOVAL TOOL		CI-7526

\* This unpainted colour cap is intended for hypo-allergenic patients who could have allergic reactions to paint. Made to order, please contact your local AB representative prior to ordering.

# SLIM HP MAGNETS & STRENGTH

SLIM HP STANDARD MAGNET 1	CI-7524-001
SLIM HP STANDARD MAGNET 2	CI-7524-002
SLIM HP STANDARD MAGNET 3	CI-7524-003
SLIM HP STANDARD MAGNET 4	CI-7524-004
SLIM HP STANDARD MAGNET 5	CI-7524-005
SLIM HP STANDARD MAGNET 6	CI-7524-006
SLIM HP STANDARD MAGNET 7	CI-7524-007

SLIM HP 3D MAGNET 1	CI-7525-001
SLIM HP 3D MAGNET 2	CI-7525-002
SLIM HP 3D MAGNET 3	CI-7525-003
SLIM HP 3D MAGNET 4	CI-7525-004
SLIM HP 3D MAGNET 5	CI-7525-005
SLIM HP 3D MAGNET 6	CI-7525-006
SLIM HP 3D MAGNET 7	CI-7525-007



If Slim HP Standard Magnet  $\geq 6$   
 ➔ Slim HP Colour Cap Large\*



MINIMUM

MAXIMUM

Increasing magnetic strength



If Slim HP 3D Magnet  $\geq 5$   
 ➔ Slim HP Colour Cap Large\*




\*The Slim HP Colour Cap Large is not indicated for use with children under the age of 3 years.

## SLIM HP MAGNET KIT

The Slim HP Magnet Kit includes a complete set of headpiece magnets, of both kinds and of every available strength. In details:

Magnet Strength	1	2	3	4	5	6	7
Standard	2 pcs	2 pcs	3 pcs	3 pcs	2 pcs	2 pcs	2 pcs
3D	2 pcs	2 pcs	3 pcs	3 pcs	2 pcs	2 pcs	2 pcs

SLIM HEADPIECE MAGNET KIT CI-7528

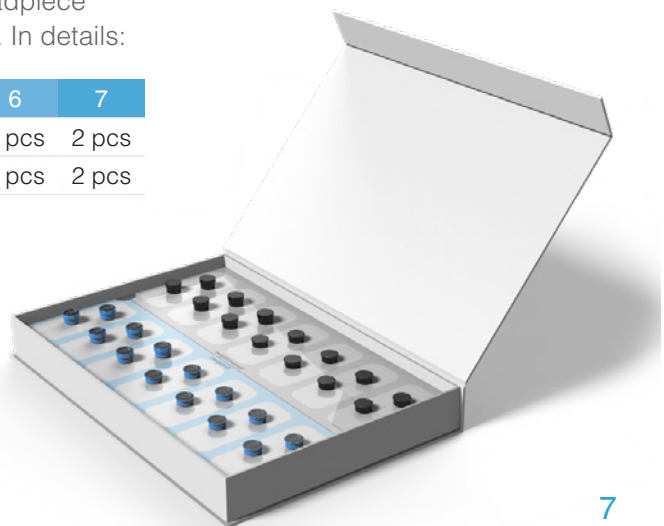



# M WATERPROOF BATTERY

The **M Waterproof Battery** gives you the same great hearing while you enjoy even the most rugged activities, all while providing secure, worry-free wearing. And AB's **AquaMic™** headpieces are the world's only waterproof sound processor microphones. This lets you enjoy an intense workout or pool party, keeping you a part of the conversation in either setting.

## COMPONENTS

1. Sound Processor
2. Cable Connector
3. Security Lock
4. Slim HP AquaMic
5. M Waterproof Battery Clip




M WATERPROOF BATTERY	CI-5068
M WATERPROOF BATTERY CLIP	CI-7323
M WATERPROOF BATTERY O-RINGS	CI-7324
M ARMBAND STANDARD (42 cm)	CI-7435-001
M ARMBAND SMALL (26 cm)	CI-7435-002



## SLIM HP AQUAMIC

SLIM HP AQUAMIC 30.5 CM (12") BEIGE	●	CI-5321-001
SLIM HP AQUAMIC 30.5 CM (12") BLACK	●	CI-5321-002
SLIM HP AQUAMIC 30.5 CM (12") BROWN	●	CI-5321-003
SLIM HP AQUAMIC 30.5 CM (12") GREY	●	CI-5321-004
SLIM HP AQUAMIC 46 CM (18") BEIGE	●	CI-5321-005
SLIM HP AQUAMIC 46 CM (18") BLACK	●	CI-5321-006
SLIM HP AQUAMIC 46 CM (18") BROWN	●	CI-5321-007
SLIM HP AQUAMIC 46 CM (18") GREY	●	CI-5321-008
SLIM HP AQUAMIC 61 CM (24") BEIGE	●	CI-5321-009
SLIM HP AQUAMIC 61 CM (24") BLACK	●	CI-5321-010
SLIM HP AQUAMIC 61 CM (24") BROWN	●	CI-5321-011
SLIM HP AQUAMIC 61 CM (24") GREY	●	CI-5321-012
SLIM HP AQUAMIC 107 CM (42") BEIGE	●	CI-5321-013
SLIM HP AQUAMIC 107 CM (42") BLACK	●	CI-5321-014
SLIM HP AQUAMIC 107 CM (42") BROWN	●	CI-5321-015
SLIM HP AQUAMIC 107 CM (42") GREY	●	CI-5321-016

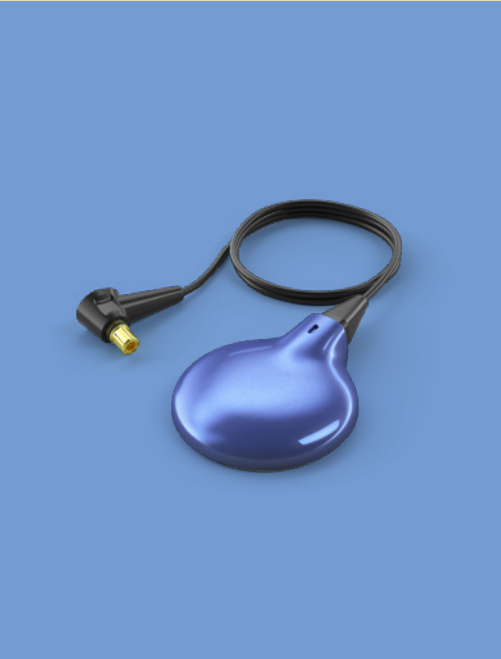

## POWER OPTIONS

You or your child will be wearing the Marvel CI sound processor nearly around the clock. This calls for a device that is comfortable, discreet, and which will not interfere with your favourite activities.

Marvel CI has a range of new battery options to suit all lifestyles, providing a **unique ergonomic housing** for a comfortable and secure fit: with the M Batteries, Small, Medium and Large, the bottom of the processor curves against the contour of the ear, keeping the device secure and conveniently tucked away. The new M Waterproof Battery is small, lightweight and ideal for swimming and outdoor activities.

The graph below illustrates the weight, the capacity and the average operating hours of each M Battery option.



Average  
Operating  
Hours

>10 hrs

>16 hrs

>16 hrs

>16 hrs

>16 hrs

>18 hrs

Weight\*

4.6 g

6.4 g

8.4 g

5.4 g

6.4 g\*\*

29.4 g

\* Processor alone: 4.6 g

\*\* With two batteries



## M BATTERY SMALL

ALPINE WHITE	●	CI-5551-110
SAND BEIGE	●	CI-5551-120
CHESTNUT	●	CI-5551-130
SILVER GREY	●	CI-5551-140
VELVET BLACK	●	CI-5551-150
CARIBBEAN PIRATE	●	CI-5551-190
PRECIOUS PINK	●	CI-5551-250
BLUE OCEAN	●	CI-5551-260
MAJESTY PURPLE	●	CI-5551-270
LAVA RED	●	CI-5551-280
RESIN BEIGE*	●	CI-5551-240

## M BATTERY LARGE

ALPINE WHITE	●	CI-5553-110
SAND BEIGE	●	CI-5553-120
CHESTNUT	●	CI-5553-130
SILVER GREY	●	CI-5553-140
VELVET BLACK	●	CI-5553-150
CARIBBEAN PIRATE	●	CI-5553-190
PRECIOUS PINK	●	CI-5553-250
BLUE OCEAN	●	CI-5553-260
MAJESTY PURPLE	●	CI-5553-270
LAVA RED	●	CI-5553-280
RESIN BEIGE*	●	CI-5553-240

## M BATTERY STANDARD

SAND BEIGE	●	CI-5555-120
CHESTNUT	●	CI-5555-130
SILVER GREY	●	CI-5555-140
VELVET BLACK	●	CI-5555-150
CARIBBEAN PIRATE	●	CI-5555-190
PRECIOUS PINK	●	CI-5555-250
BLUE OCEAN	●	CI-5555-260
MAJESTY PURPLE	●	CI-5555-270
RESIN BEIGE*	●	CI-5555-240

## M BATTERY MEDIUM

ALPINE WHITE	●	CI-5552-110
SAND BEIGE	●	CI-5552-120
CHESTNUT	●	CI-5552-130
SILVER GREY	●	CI-5552-140
VELVET BLACK	●	CI-5552-150
CARIBBEAN PIRATE	●	CI-5552-190
PRECIOUS PINK	●	CI-5552-250
BLUE OCEAN	●	CI-5552-260
MAJESTY PURPLE	●	CI-5552-270
LAVA RED	●	CI-5552-280
RESIN BEIGE*	●	CI-5552-240

## M ZN-AIR BATTERY PAK

ALPINE WHITE	●	CI-5501-110
SAND BEIGE	●	CI-5501-120
CHESTNUT	●	CI-5501-130
SILVER GREY	●	CI-5501-140
VELVET BLACK	●	CI-5501-150
CARIBBEAN PIRATE	●	CI-5501-190
PRECIOUS PINK	●	CI-5501-250
BLUE OCEAN	●	CI-5501-260
MAJESTY PURPLE	●	CI-5501-270
LAVA RED	●	CI-5501-280
RESIN BEIGE*	●	CI-5501-240

## CHARGER

M BATTERY CHARGER	CI-5607
M BATTERY WALLET	CI-7436



\*These unpainted batteries are intended for hypo-allergenic patients who could have allergic reactions to paint. Made to order, please contact your local AB representative prior to ordering.



## M EARHOOK

M EARHOOK SMALL	CI-5711-100
M EARHOOK STANDARD	CI-5711-110
M EARHOOK WITH GRIP SMALL	CI-5711-200
M EARHOOK WITH GRIP STANDARD	CI-5711-210
M PIN REMOVAL TOOL (2-pack)	CI-6147-014
M EARHOOK PIN (10-pack)	CI-6147-015



## M ACOUSTIC EARHOOK

M ACOUSTIC EARHOOK RIGHT 00	CI-5851-001
M ACOUSTIC EARHOOK RIGHT 0	CI-5851-002
M ACOUSTIC EARHOOK RIGHT 1	CI-5851-003
M ACOUSTIC EARHOOK RIGHT 2	CI-5851-004
M ACOUSTIC EARHOOK RIGHT 3	CI-5851-005
M ACOUSTIC EARHOOK LEFT 00	CI-5851-006
M ACOUSTIC EARHOOK LEFT 0	CI-5851-007
M ACOUSTIC EARHOOK LEFT 1	CI-5851-008
M ACOUSTIC EARHOOK LEFT 2	CI-5851-009
M ACOUSTIC EARHOOK LEFT 3	CI-5851-010
M ACOUSTIC EARHOOK TOOL KIT	CI-6150



## M T-MIC MICROPHONE

M T-MIC SMALL	CI-5837-100
M T-MIC MEDIUM	CI-5837-200
M T-MIC LARGE	CI-5837-300
M T-MIC PINS (25 pieces)	CI-6147-005
M T-MIC COVERS (25 pieces)	CI-6147-006
M T-MIC COVERS (8 pieces)	CI-6147-007
M PIN REMOVAL TOOL (2-pack)	CI-6147-014

## WEARING & RETENTION OPTIONS

Naída CI M and Sky CI M feature a variety of wearing options that can be customised to the changing needs of you and your child.

At Advanced Bionics, we know little ears best. That's because we've spent years working to ensure our technology stays where it should be—out of your child's way as they explore the world.

M WEARING CLIP	CI-7322
M SNUGGIE	CI-7325
M RETENTION CUFF SMALL (44 mm)	CI-7326-001
M RETENTION CUFF LARGE (64 mm)	CI-7326-002

Nothing on the ear  
M WEARING CLIP



Secure on the ear  
M SNUGGIE

Secure behind the ear  
M RETENTION CUFF





# PHONAK SOUND DELIVERY SYSTEM 4.0

VENTED DOME 4.0 S (package)	054-0809
VENTED DOME 4.0 M (package)	054-0810
VENTED DOME 4.0 L (package)	054-0811
POWER DOME 4.0 S (package)	054-0820
POWER DOME 4.0 M (package)	054-0821
POWER DOME 4.0 L (package)	054-0822
CERUSHIELD™ DISK	098-0445
MEASUREMENT TOOL FOR RECEIVER CABLE LENGTH	026-0697



Power Dome 4.0



Vented Dome 4.0



CeruShield Disk

## ADDITIONAL ACCESSORIES

### Programmaming

NOAHLINK WIRELESS PROGRAMMING INTERFACE	058-0139*
M PROGRAMMING CABLE	CI-5826
CPI-3 USB CABLE M SERIES (2M)	CI-6408-002

### Device Check

M LISTENING CHECK	CI-5827
LISTENING CHECK EARBUDS	CI-5822

### Kits & Cases

DRYING KIT	CI-7313-001
M CARRYING CASE	305-Q144*

\*SAP number only

# AB REMOTE APP FOR BETTER CONTROL

Being able to easily control your sound processor means you or your child will always get the best sound possible. Our AB Remote app makes this easy and discreet by providing convenient management of all your hearing devices on one smartphone. The AB Remote app is available for free on both iOS and Android devices.







## BUILT-IN CONNECTIVITY

Unlike any other cochlear implant platform, Naída CI M and Sky CI M offer the ability to stream calls and audio directly from virtually every smartphone and Bluetooth®-enabled device. Our **Phonak RogerDirect™ technology** boosts hearing performance in challenging listening situations, without extra devices attached to the sound processor.

**RogerDirect** is a Marvel feature that allows the Roger microphones to stream directly to the hearing device without having to attach an external receiver. Compared to AB's previous generation, the Marvel CI sound processors provide full Roger Performance & Connectivity being **23% shorter** and **34% lighter**.



## ROGER ACCESSORIES

ROGER X (02)		052-3113-X02P5
ROGER INSTALLER		055-3033
ROGER SELECT CHAMPAGNE		052-3471-000P5
ROGER SELECT GRAPHITE GREY		052-3471-000P7
ROGER SELECT PEARL WHITE		052-3471-000W3
ROGER SELECT <b>iN</b> CHAMPAGNE		056-3002-P5
ROGER SELECT <b>iN</b> GRAPHITE GREY		056-3002-P7
ROGER SELECT <b>iN</b> PEARL WHITE		056-3002-W3

# PHONAK ACCESSORIES

PHONAK TV CONNECTOR D <sup>1</sup> (2.0)	076-3006-0612
PHONAK REMOTECONTROL	076-0065-P5
PHONAK PARTNERMIC	076-4001-P510

1. As long as stock lasts



ROGER X



ROGER INSTALLER



ROGER SELECT  
ROGER SELECT iN



PHONAK  
TV CONNECTOR



PHONAK  
REMOTECONTROL



PHONAK  
PARTNERMIC

## BIMODAL SOLUTION

At Advanced Bionics, we're working every day to ensure that your and your child's unique hearing needs are met. Together with Phonak, we're part of the largest and most innovative hearing company dedicated to delivering solutions for patients with every type and degree of hearing loss.

If you use a cochlear implant and a hearing aid, you can have a Naída CI M or a Sky CI M on one ear and a **Phonak Naída™ Link M** or a **Phonak Sky™ Link M** hearing aid on the other ear. This customised approach is designed to ensure both devices will process sounds in the same way. Together, they provide consistent sound experiences and excellent speech understanding, unlike having a hearing aid and a cochlear implant that are not designed to work together.







## Designed to work together

For individuals with some hearing in one ear, wearing the right hearing aid helps them to hear even better than with a CI alone. Wearing a Naída Link M together with Naída CI M gives you a consistent sound impression between ears, and allows you to adjust both devices with a single control.

The same is true for Sky Link M and Sky CI M.

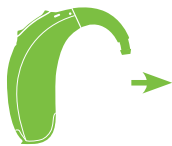



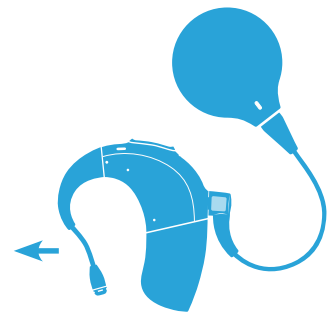





Naída Link M / Sky Link M  
hearing aid



Naída CI M / Sky CI M  
sound processor



## PHONAK NAÍDA LINK M & PHONAK SKY LINK M

PHONAK NAÍDA LINK M SAND BEIGE		050-0801-P1
PHONAK NAÍDA LINK M CHESTNUT		050-0801-P4
PHONAK NAÍDA LINK M SILVER GREY		050-0801-P6
PHONAK NAÍDA LINK M VELVET BLACK		050-0801-P8
PHONAK SKY LINK M SAND BEIGE		050-0802-P1
PHONAK SKY LINK M CHESTNUT		050-0802-P4
PHONAK SKY LINK M VELVET BLACK		050-0802-P8
PHONAK SKY LINK M CARIBBEAN PIRATE		050-0802-Q3
PHONAK SKY LINK M PRECIOUS PINK		050-0802-T3



# PROGRAMMEMES AND FEATURES OVERVIEW

PROGRAMME / FEATURE		NAÍDA CI M90	SKY CI M90
<b>AUTOSENSE OS 3.0 / AUTOSENSE SKY OS 3.0</b>			
Calm Situation	Use for speech understanding and listening comfort in 1:1 conversations and for quiet environments.	●	●
Speech in Noise	Use for speech understanding and comfort when in noisy environments and groups.	●	●
Comfort in Noise	Use for listening comfort in noisy situations when no speech is present.	●	●
Music	Programme settings for a more natural music experience.	●	●
Speech in Loud Noise	Use for speech understanding when in noisy and crowded environments. The StereoZoom feature in the programme allows a particular focus on the sound source in front of the listener.	●	●
Speech in Car	Optimised settings designed for a stable listening environment in the car (awareness).	●	●
Comfort in Echo	Use this programme for comfortable speech understanding in environments that have a lot of echo/reverberation.	●	●
Streamed Speech	Used when streaming speech from A2DP or AirStream.	●	●
Streamed Music	Used when streaming music from A2DP or Airstream.	●	●
<b>ADDITIONAL MANUAL PROGRAMMES (max 4)</b>			
Speech in 360°	This programme enables the selection of the listening direction in situations where the sound source cannot be faced. Focus direction can be fixed or automatic.	●	●
DuoPhone	The key feature DuoPhone (M90/Sky only) streams the signal in real-time to the other ear so the caller's voice is heard in both ears without the need for any additional accessories.	●	●
Off-Ear Programme	Use this programme when the device is configured for off-ear wearing using a headpiece microphone.	●	●
Custom Programme	–	●	●
<b>STREAMING PROGRAMMES</b>			
RogerDirect + mic	Combines input from Roger microphone source and device microphones.	●	●
Phone call + mic	Combines input from streaming phone source and device microphones.	●	●

PROGRAMME / FEATURE		NAÍDA CI M90	SKY CI M90
<b>SIGNATURE FEATURES</b>			
UltraZoom	UltraZoom is the adaptive multi-channel beamforming technology, which improves speech intelligibility within noisy situations.	●	●
UZ + SNR-Boost	Designed to improve SNR in an UltraZoom programme by providing further attenuation of sounds from the back.	●	●
StereoZoom	StereoZoom is an adaptive, multi-channel beamforming technology that adapts to moving and asymmetric noise sources. It requires compatible bilateral devices.	●	●
DuoPhone	Automatically streams the caller's voice to the other ear, so that the voice is heard in both ears simultaneously.	●	●
Real Ear Sound	Real Ear Sound simulates the performance characteristics of the T-Mic.	●	●
WindBlock	WindBlock utilizes a wind identification and classification system. The powerful, adaptive algorithm automatically suppresses wind noise.	●	●
EchoBlock	EchoBlock is designed for reverberant environments like large assembly halls, places of worship, and rooms with reflective surfaces. EchoBlock works by instantly removing the disturbing reflections of the reverberant sound.	●	●
SoundRelax	SoundRelax instantly and accurately identifies impulse sounds and applies gain reduction relative to the original peak signal characteristics. SoundRelax cushions otherwise unpleasant impulse sounds.	●	●
QuickSync	One-touch controls two. The instantaneous broadband transfer of data between two instruments ensures that volume and programme settings are always balanced and synchronised. One-touch on the programme selection button or volume control results in an instant identical change in the other instrument.	●	●
Environmental balance	Volume modification capability that allows clients to adjust the volume of the streamed signal versus the surrounding environmental sounds by the press of the volume control.	●	●



## **ADVANCED BIONICS LLC**

28515 Westinghouse Place  
Valencia, CA 91355, United States

T: +1.877.829.0026

T: +1.661.362.1400

F: +1.661.362.1500

*info.us@advancedbionics.com*

## **ADVANCED BIONICS AG**

Laubisrütistrasse 28  
8712 Stäfa, Switzerland

T: +41.58.928.78.00

F: +41.58.928.78.90

*info.switzerland@advancedbionics.com*

## **ADVANCED BIONICS UK & EIRE**

2 Breaks House, Mill Court, Great Shelford,  
Cambridge, CB22 5LD, UK

T: +44.1223.847.888

F: +44.1223.847.898

*info.uk@advancedbionics.com*

For information on additional AB locations, please visit  
[advancedbionics.com/contact](https://www.advancedbionics.com/contact)

Advanced Bionics – A Sonova brand

Please contact your local AB representative for regulatory approval and availability in your region.

\*Not all colours and sizes are available in all markets.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Sonova AG is under license.

Apple, the Apple logo, iPhone, and iOS are trademarks for Apple Inc., registered in the US and other countries. App Store is a service mark of Apple Inc.

Android, Google Play, and the Google Play logo are trademarks of Google Inc.